

## A-FLX™ Animal Tissue FFPE Wistar Rat Heart

### GENERAL INFORMATION

Product Name: A-FLX™ Animal Tissue  
FFPE Wistar Rat Heart

Reference Number: 7020-0410 Block  
7020-0420 Slide (4µm)

Date of Manufacturing: See product label

Lot Number: See product label

Intended Use: For research use only

### DESCRIPTION

Species: Rat

Strain: Wistar

Tissue: Heart

Fixation Condition: 10% neutral buffered formalin (NBF) for 24 hours at 24-27°C

Product Format:

Block: Paraffin embedded block

Slide: One unstained section mounted on Superfrost™ Plus slide. Section thickness: 4µm

### IMAGE

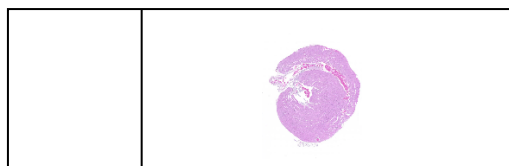


Illustration of an FFPE slide with  
Hematoxylin & Eosin (H&E) staining image

### STORAGE CONDITION

2-8 ° C with desiccation (Recommended)

### \*RECOMMENDED STABILITY RE-TESTING SCHEDULE:

Block	5 years
Slide	1 year

\*FFPE sample stability is not universally established and is biomarker and assay dependent. More frequent re-testing may be needed for certain labile biomolecules.

### RECOMMENDED PROCEDURES

#### Staining using FFPE slides:

1. Bake slides at 60°C for 30-60 min.
2. Deparaffinize 3 times in Xylene or Xylene substitute for 5 min each time.
3. Rinse 2 times in 100% ethanol for 1 min each.
4. Rehydrate in ethanol series (95% 1 min, 70% 1 min, distilled H<sub>2</sub>O 2 times for 1 min each).
5. Proceed to staining protocol.

### SAFETY AND PRECAUTIONS

This product does not contain hazardous material. Wear appropriate personal protective equipment (PPE) when handling reagents and biological specimens.



#### Acepix Biosciences, Inc.

32990 Alvarado-Niles Road, Suite 910  
Union City · CA · 94587 · USA

Tel: 1 (888) 866-0580 · Web: [www.acepibio.com](http://www.acepibio.com)